



# Roberts Lab\_2011-02-17 11-30-30\_CC009827.pcrd

2/17/2011 2:47 PM

## Report Information

User: BioRad\Roberts Lab  
Data File Name: Roberts Lab\_2011-02-17 11-30-30\_CC009827.pcrd  
Data File Path: C:\Users\srlab\Documents\My Dropbox\Roberts Lab CFX96 Data (7)\Sam  
Selected Well Group: All Wells

## Experiment Setup

### Run Information

Run Date: 2/17/2011 11:30:39 AM  
Run User: BioRad\Roberts Lab  
ID:  
Notes:  
Sample Volume: 25  
Temperature Control Mode: Calculated  
Lid Temperature: 105  
Base Serial Number: CC009827  
Optical Head Serial Number: 785BR3659

### Protocol

1: 95.0°C for 10:00  
2: 95.0°C for 0:15  
3: 55.0°C for 0:15  
4: 72.0°C for 0:15  
Plate Read  
5: GOTO 2, 39 more times  
6: 95.0°C for 0:10  
7: Melt Curve 65°C to 95°C : Increment 0.5°C for 0:05  
Plate Read


### Plate Display

	1	2	3	4	5	6	7	8	9	10	11	12
A	Link 18s BB01											
B	Link 18s BB09											
C	Pos 18s BB16.gDNA											
D	NTC 18s											
E	NTC 18s											
F												
G												
H												

### Quantitation

Step #: 4  
Analysis Mode: Baseline Subtracted Curve Fit  
Ct Determination: Single Threshold  
Baseline Method per Fluorophore:  
FAM: Auto Calculated  
Threshold Setting per Fluorophore:  
FAM: 474.65, Auto Calculated

---

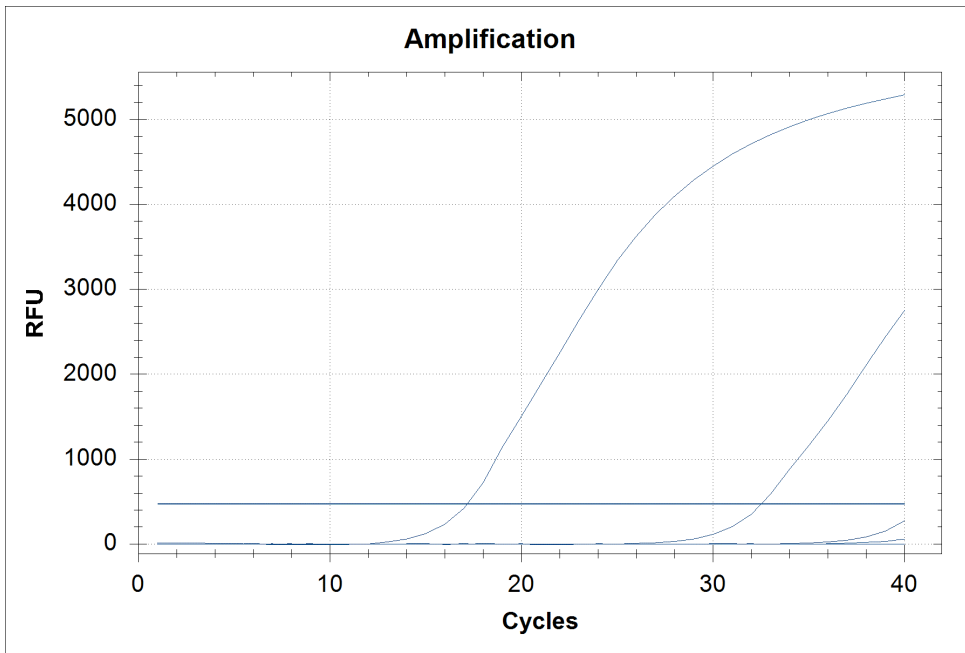



#### Quantitation Data

Well	Fluor	Content	Target	Sample	Threshold Cycle ( C(t) )	C(t) Mean	C(t) Std. Dev
A01	FAM	Unkn	18s	BB01	32.52	32.52	0.000
B01	FAM	Unkn	18s	BB09	N/A	0.00	0.000
C01	FAM	Pos Ctrl	18s	BB16 gDNA	17.18	17.18	0.000
D01	FAM	NTC	18s		N/A	0.00	0.000
E01	FAM	NTC	18s		N/A	0.00	0.000

#### Melt Curve

Step #: 7

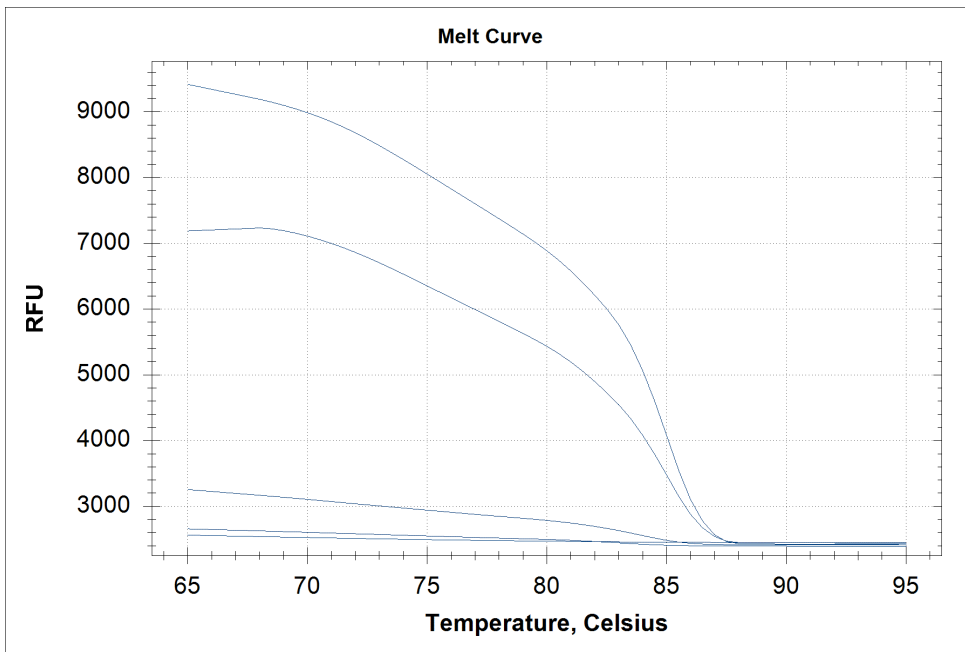

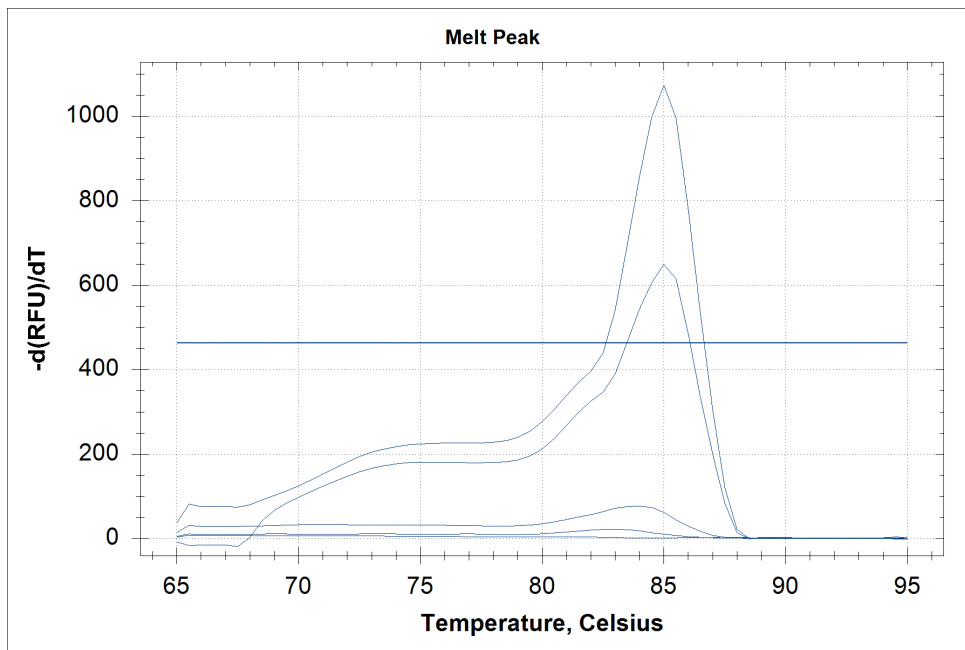

**Melt Curve Data**

Well	Fluor	Content	Sample	Melt Temp
A01	FAM	Unkn	BB01	85.00
C01	FAM	Pos Ctrl	BB16 gDNA	85.00